




# Mirage Kidsta™ Nasal Mask

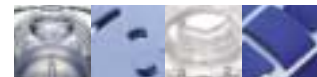

**The reassuring mask of choice for sleep-disordered breathing (SDB) in childhood**

ResMed is the world leader in mask technology for SDB and we are bringing that expertise to a mask designed for children seven years and older.

The Mirage Kidsta™ Nasal Mask is based on the revolutionary range of adult Mirage™ designs. It combines a number of proven, practical and popular benefits with some child-friendly modifications. This ResMed mask is designed specifically for children who need continuous positive airway pressure (CPAP) or bilevel therapy.

ResMed's Kidsta Nasal Mask—with the proven and trusted Mirage dual-wall cushion technology—is a premium mask that combines comfort and performance in an unobtrusive design.

**Stabilizing panel**  
holds headgear in correct position for optimal seal

**Headgear**  
provides set and forget convenience

**Mask elbow**  
connects tubing to mask frame. Provides greater ease for most children to click-on and click-off

**Angle selector**  
enables customized fit for optimal comfort

**Easy locking mechanism**  
provides easier handling and enables headgear clip to easily lock into mask frame



## Technology highlight—leaders in mask technology

ResMed is renowned for its pioneering technologies and commitment to producing world-class products.

The Mirage Kidsta Nasal Mask was designed and manufactured based on anthropometric (facial and head) measurements on a range of children.

Data was assimilated and combined with ResMed's unique Mirage dual-wall cushion technology to help develop a mask cushion and headgear specifically designed for a child.



## Benefits

- Dual-wall cushion reduces pressure on small noses and optimizes seal and comfort.
- Easy click-on, click-off parts are ideal for children's hands to grasp.
- Unique instruction booklet, specifically designed to support parent and child.
- Angle selector and stabilizing panels provide customized fit, enhanced stability, seal and comfort.
- Unobtrusive, streamlined design allows eye contact and clear vision. The open design also reduces feelings of anxiety and claustrophobia.
- Convenient set and forget headgear minimizes the need for continuous adjustment and allows greater control for smaller hands.
- Child can wear mask without elbow attached to tubing to facilitate acclimatization to therapy.



Open design provides child with a clear field of vision and encourages acclimatization.



Child-friendly design allows smaller hands to grasp parts easily.



The Mirage Kidsta Nasal Mask comes with a child's instruction booklet introduced by the mask mascot Oke. The guide includes a diary to help chart progress and is valuable when it comes to acclimatization and compliance.

### Mirage Kidsta Nasal Mask System

**Product Code** 61011